


# BiPAP A40™ Oxygen Therapy



**Pressure Valve**  
Respironics Part  
302418

**Oxygen Adapter**  
AirLife Cat. 004081

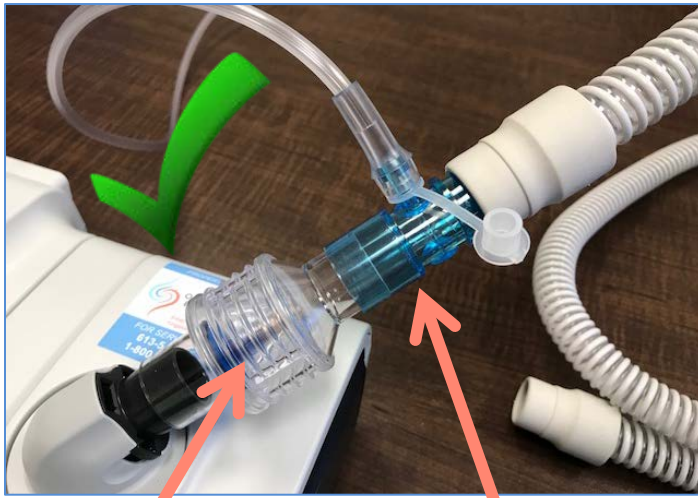

Oxygen therapy may be used with the BiPAP A40™.

To add oxygen, you will need:

- Pressure Valve
- Oxygen Adapter

The pressure valve directs oxygen flow to the patient. It is important to stop oxygen from flowing back into the machine.<sup>1</sup>

The oxygen adapter has a connection for low flow oxygen tubing.

Once the pressure valve and oxygen adapter are connected as shown on the left, attach the oxygen tubing from your oxygen supply system and set flow as prescribed (L/min).

Oxygen Flow  $\leq 10$  L/min<sup>1</sup>

Do Not Connect  
the Oxygen  
Adapter to the  
machine outlet!



## References

1. Philips Respironics™ BiPAP A40™ User Manual

## Disclaimer

This information is provided to you for educational purposes only and should not at any time be considered for a substitute for professional advice from your physician or other qualified health care professional. The VEP does not assume any responsibility for application of the content of this document or for any information that may be present in the documents or websites cited as a reference. These references are provided for informational purposes only and do not represent VEP endorsement of any companies or products. Always consult your physician and/or a qualified health care professional to learn more about recommendations specific to your health needs.



640 Cataraqui Woods Drive | Unit 6 | Kingston, Ontario | K7P 2Y5  
Tel: 613-548-6156 | Toll-Free in Ontario: 1-800-633-8977 |  
Fax: 613-548-6157  
Website: [www.ontvep.ca](http://www.ontvep.ca)