

## How to Download Data from the USB Stick and Send by Email:

### Stellar 150™




Insert an **empty** USB stick into either USB port found at the back of the Stellar 150™. Only 1 USB stick can be used at a time.



After inserting the USB stick the following window comes up.

Using the scroll wheel, select **write Set/Log** and click the scroll wheel to Start.

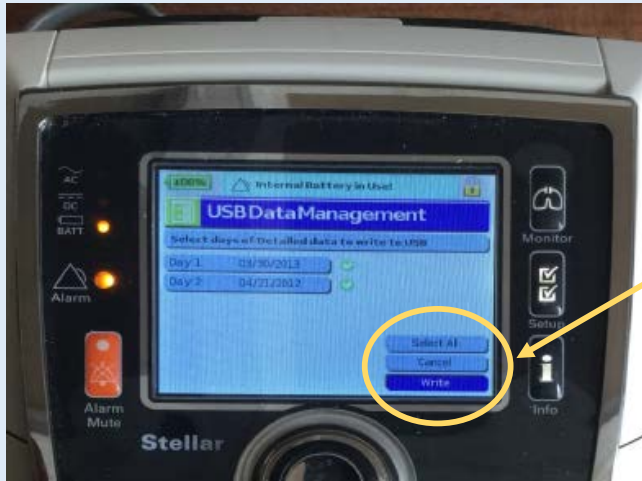


The following window will come up:

Using the scroll wheel, click on **Select All** for the days of detailed data, and click **write**.



Wait for the progress bar to complete (this may take several minutes)

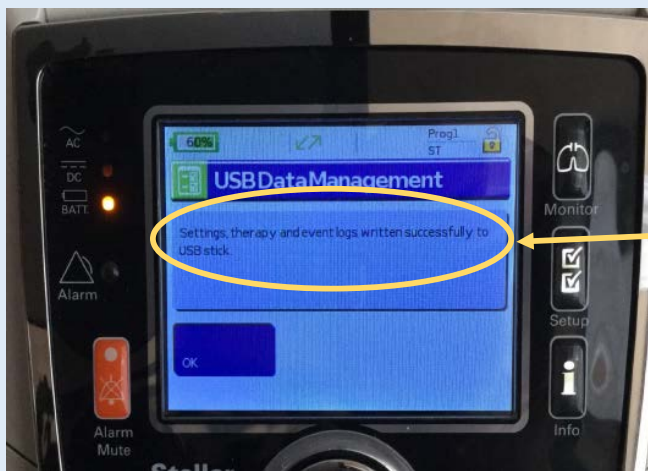


When you see the message "**Settings, therapy and event logs written successfully to USB stick**" the data has been successfully copied to the USB stick. Push the scroll wheel to hit "OK"


## To Download and Send the File



Take the USB stick that you used to download the data from the Stellar, and insert it into USB port of your PC.

Create a folder on your PC (**make sure you remember the location of the folder**), and name it: **First Name, Last Name & date of download**

 John Doe April 20, 2020

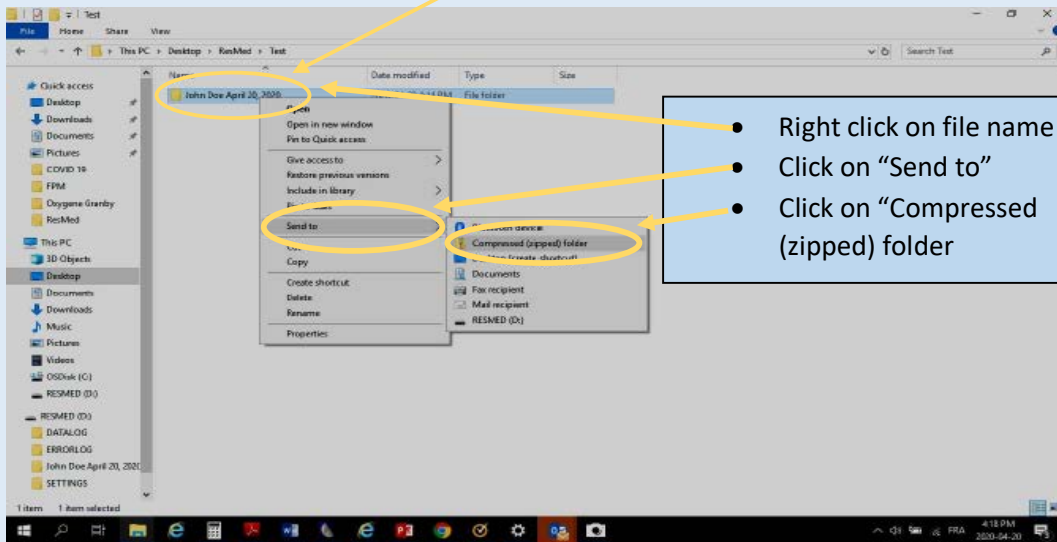

- Right click on file name
- Click on “Send to”
- Click on “Compressed (zipped) folder”

- Copy all the files and folders from the USB stick where you did the download, to the folder you just created on your PC.
- **Zip that folder and send the zipped folder by email to destination.**
- Email to: **vep@kingstonhsc.ca**

## **To Receive and Download the Emailed File**

- Open up the email with the zipped folder and unzip the patient's folder.
- Copy all the files and folders within that unzipped folder onto an empty USB stick.
- Open Rescan™ software and download the data from the USB stick.
- You will need to create a new patient in ResScan™ or add the data to an existing patient.

### **Disclaimer**

This information is provided to you for educational purposes only and should not at any time be considered for a substitute for professional advice from your physician or other qualified health care professional. The VEP does not assume any responsibility for application of the content of this document or for any information that may be present in the documents or websites cited as a reference. These references are provided for informational purposes only and do not represent VEP endorsement of any companies or products. Always consult your physician and/or a qualified health care professional to learn more about recommendations specific to your health needs.



640 Cataraqui Woods Drive, Unit 6, Kingston, Ontario, K7P 2Y5  
Tel: 613-548-6156, Toll-Free in Ontario: 1-800-633-8977, Fax: 613-548-6157  
Web: [www.ontvep.ca](http://www.ontvep.ca)